

HITACHI OASIS

LEGEND

- RF SHIELD REQUIRED. ARCHITECT TO VERIFY WALL CONSTRUCTION WITH RF VENDOR. IT MUST PROVIDE 80 DB ATTENUATION OR BETTER FROM 1-70MHZ. (TYPICAL 5" RF WALL THICKNESS REFLECTED.)
- RF VIEW WINDOW. TYPICALLY 4'W X 3'H. (TYP. 36" AFF)
- REMOVABLE EQUIPMENT ACCESS PANEL FOR DELIVERY OF THE SYSTEM. MINIMUM 8'-6"W X 8'-5"H.
- JUNCTION BOX FOR POWER TO RFP CABINET.
- MANUAL DISCONNECT
- 480V, 30 AMP, 3PH, 4 WIRE. HUBBELL L16-30 RECEPTACLE (OR EQUAL) FOR RAMP-UP TOOL USAGE.
- 115V, 15 AMP CONVENIENCE OUTLETS. MINIMUM OF TWO (2) PER ROOM REQUIRED FOR SERVICE AND MAINTENANCE.
- PHONE LINES - ONE (1) REQUIRED IN EQUIPMENT ROOM AND ONE (1) REQUIRED IN CONTROL ROOM
- DATA PORT FOR BROAD BAND NETWORK CONNECTION TO STANDARD 16 BASE T LINE. MAY BE THROUGH CUSTOMERS IN-HOUSE NETWORK.
- 6" X 24" ON-WALL WIREWAY WITH REMOVABLE COVER AND THREE (3) DIVIDERS. METAL ONLY. MUST BE NON-FERROUS IN THE SCAN ROOM.
- 6" X 24" ON-FLOOR WIREWAY WITH REMOVABLE COVER AND THREE (3) DIVIDERS. NON-FERROUS METAL ONLY.
- 3" X 6" OVERHEAD CABLE TRAY (WALKER SPMA-A-6-3S OR EQ). NON-FERROUS METAL ONLY.
- MAGNET ISOCENTER

EQUIPMENT BY HMSA


- (A) OASIS GANTRY
- (B) PATIENT TABLE
- (C) OPERATOR'S WORKSTATION
- (D) RFP CABINET
- (E) GCFA CABINET
- (F) MCU FILTER PANEL
- (G) EMERGENCY RUN-DOWN BUTTON (ERDU)
- (H) RFP FILTER PANEL
- (J) COIL STORAGE CABINET
- (K) HELIUM COMPRESSOR
- (M) SENSE UNIT
- (N) HEAT EXCHANGER
- (P) R1200 OUTDOOR CHILLER

BY OTHERS

- (2) 6" X 24" NON-FERROUS ON-WALL WIREWAY
- (3) 6" X 24" NON-FERROUS ON-FLOOR WIREWAY
- (4) A.C. POWER PANEL (SIZED PER CODE)
- (5) OXYGEN PROBE
- (6) OXYGEN MONITOR
- (7) LASER IMAGER
- (8) CHAIR
- (9) 1 3/4" WAVE GUIDES
- (10) DELIVERY ACCESS PANEL
- (11) FLOOR DRAIN
- (12) CASEWORK
- (13) HOSE BIBB
- (14) HELIUM TANK RESTRAINT
- (15) 3" X 6" OVERHEAD CABLE TRAY
- (16) 3 1/2" X 6" WIREWAY

NOTES

1. 480V, 30A, 3PH, 4 WIRE. NEMA L16-30 RECEPTACLE FOR USE WITH RAMP-UP TOOL. TO INITIALIZE MAGNETIC FIELD. THE OASIS SYSTEM DRAWS ONLY MINIMAL POWER DURING INITIALIZATION OF FIELD. RAMP-UP TOOL LOAD SHOULD NOT BE INCLUDED WHEN CALCULATING TOTAL LOAD ON PANEL FROM WHICH THE OASIS MRI SYSTEM DRAWS ITS POWER.
2. FILTER PANELS MUST BE PLACED AS SHOWN ON THIS DRAWING TO ENSURE PROPER FUNCTION OF THE OASIS MRI SYSTEM.
3. WIREWAYS MUST BE RUN AS SHOWN ON THIS DRAWING TO ENSURE PROPER FUNCTION OF THE OASIS MRI SYSTEM.
4. CUSTOMER'S ARCHITECT/CONTRACTOR IS RESPONSIBLE FOR VERIFYING ALL POWER AND HVAC SPECIFICATIONS FOR EQUIPMENT NOT SUPPLIED BY HMSA.
5. THIS DRAWING IS NOT FOR CONSTRUCTION PURPOSES. CUSTOMER IS RESPONSIBLE FOR HIRING QUALIFIED DESIGN PROFESSIONALS TO PREPARE CONSTRUCTION DOCUMENTS. CONSTRUCTION SHALL BE PERFORMED IN ACCORDANCE WITH THIS SITE SPECIFIC DRAWING. THE SPECIFICATIONS IN THE HMSA STANDARD DETAILS FOR THIS SYSTEM AND ALL APPLICABLE CODES AND REGULATIONS.
6. THE MAXIMUM WALL THICKNESS FROM THE SCAN ROOM SIDE OF THE RF SHIELD TO THE FINISHED INSIDE FACE OF THE EQUIPMENT ROOM WALL IS 11" IN ORDER TO ENSURE THE PROPER FITMENT OF THE RFP FILTER PANEL.
7. HMSA RESERVES THE RIGHT TO WITHHOLD THE APPROVAL OF A SITE PENDING AN ON-SITE REVIEW AND APPROPRIATE QUALIFICATION TESTING BY HMSA PERSONNEL.
8. CONSTRUCTION DRAWINGS ARE TO BE SUBMITTED TO THE HMSA SITE PLANNING DEPARTMENT FOR REVIEW PRIOR TO CONSTRUCTION.
9. CONTACT THE HMSA SITE PLANNING DEPARTMENT AT 800-800-3106 IMMEDIATELY REGARDING ANY DEVIATION BETWEEN THIS DRAWING, HMSA OASIS STANDARD DETAILS AND/OR CONSTRUCTION DOCUMENTS.
10. THIS DRAWING IS THE PROPERTY OF HITACHI MEDICAL SYSTEMS AMERICA, INC. (HMSA) AND IS NOT FOR CONSTRUCTION PURPOSES. THIS DRAWING, IN CONJUNCTION WITH THE STANDARD DETAILS FOR THE OASIS MRI SYSTEM, IS FOR THE USE OF THE CUSTOMER'S DESIGN PROFESSIONALS TO PREPARE APPROPRIATE CONSTRUCTION DOCUMENTS IN ACCORDANCE WITH ALL APPLICABLE CODES AND REGULATIONS. HMSA RESERVES THE RIGHT TO REUSE TO INSTALL ANY OR ALL EQUIPMENT WHEN THE ROOM PREPARATION DOES NOT MEET SPECIFICATIONS. THE CURRENT STANDARD DETAILS MAY BE OBTAINED BY CONTACTING HMSA'S SITE PLANNING DEPARTMENT AT 800-800-3106.

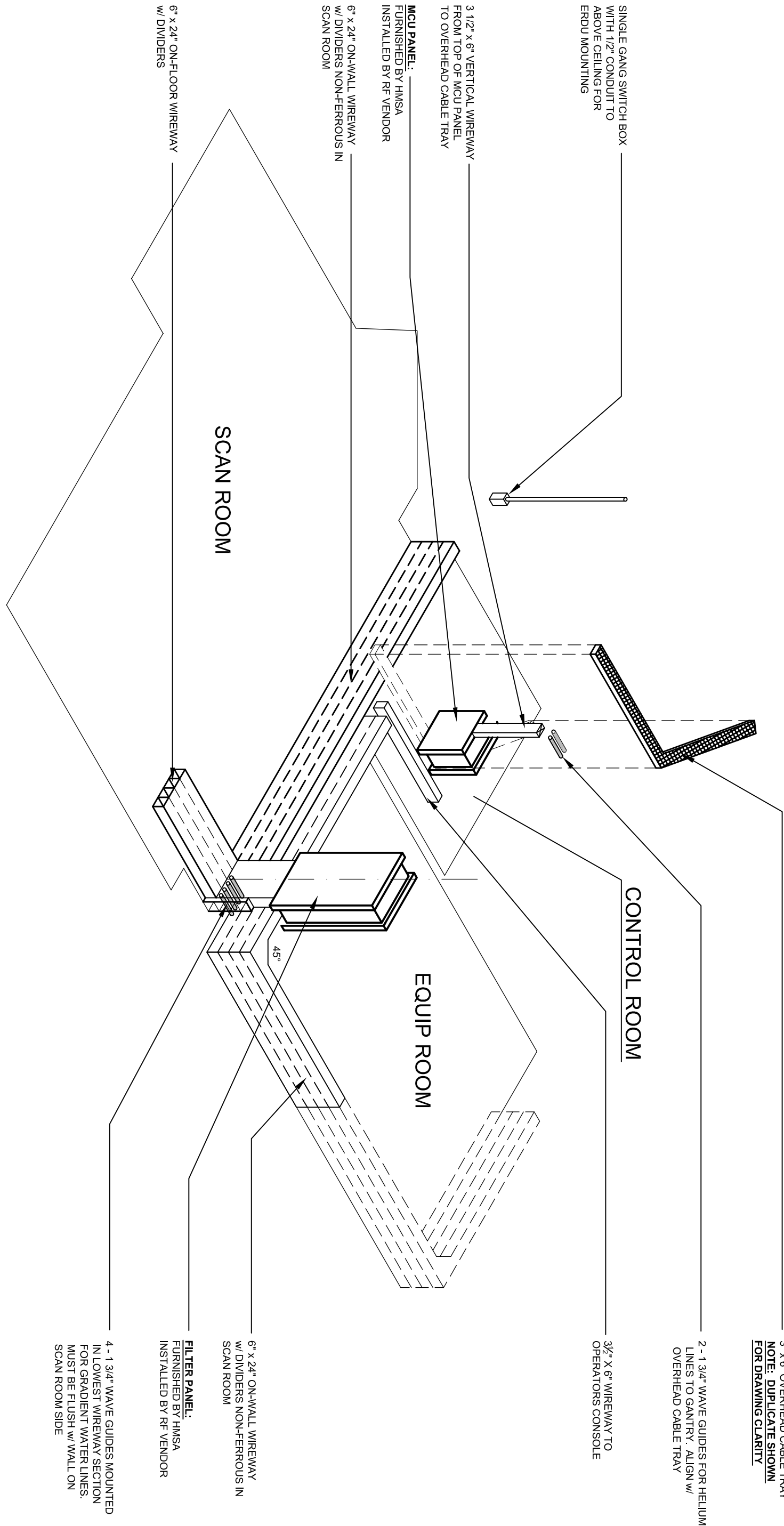

WIREWAY ISOMETRIC


SCALE: NONE

PRELIMINARY  
NOT FOR CONSTRUCTION

CONTRACTOR TO  
VERIFY DIMENSIONS &  
CONDITIONS AT THE BUILDING

CHILLER PLAN

SCALE: NONE

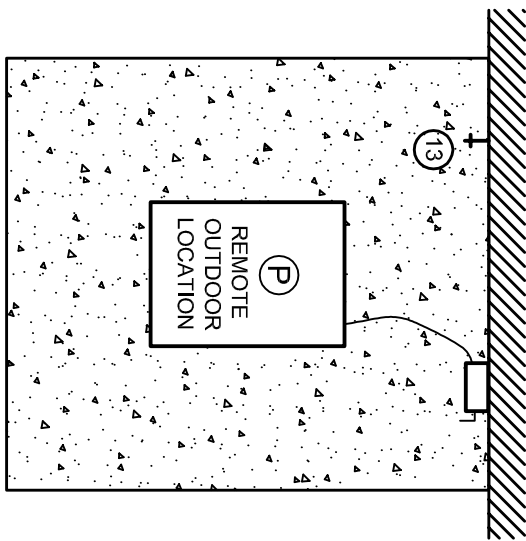

HITACHI  
Inspire the Next

HITACHI MEDICAL SYSTEMS, AMERICA  
1959 SUMMIT COMMERCE PARK TWINSBURG, OH 44087-2371  
(800) 800-3106 WWW.HITACHIMED.COM

OASIS SITE PLAN FOR:  
OVERTON BROOKS VAMC  
SHREVEPORT, LA

DATE: 08-04-2014	SCALE: AS NOTED
DRAWN BY: TJP	CHECKED BY:
FILE: 100xxx_Overton Brooks VAMC_Shreveport_LA_OASIS.dwg	
DRAWING NUMBER	: 100xxx 1 of 2