

SHIP TO:  
VA MEDICAL CENTER  
2615 E CLINTON AVE  
FRESNO, CA 93703-2223

PO#: 570-B70026

## Luminos Agile Max

**All items listed below are included for this system:** *(See Detailed Technical Specifications at end of Proposal.)*

Qty	Part No.	Item Description
1	14436607	<b>Luminos Agile Max</b> Luminos Agile Max is a highly versatile table side controlled X-ray diagnostic system meeting all requirements with respect to ease of use, image quality, connectivity and low dose. It features full tableside control of system movements and imaging parameters and allows you to stay with the patient throughout the examination such as interventional procedures (e.g. ERCP) or pediatric patients. The lifting table and the widest table/detector distance in the market are specially designed to accommodate the growing population of obese and bariatric patients. The large 43 cm x 43 cm MAX dynamic flat detector ensures excellent coverage and provides high-resolution images for more accuracy and efficiency. The syngo FLC digital one-stop workflow from patient registration to image documentation offers fast and easy operation for consistent exam settings and a consistent image impression.
1	14417363	<b>Upgrade to 80 kW PLF</b> Generator upgrade to 80kW for increased performance to extend the range of applications
1	14418476	<b>System covers</b> Cover for systems without compression device.
1	14428860	<b>Transparent grid 13/92, F115</b> Highly selective anti-scatter grid for scattered radiation reduction.
1	04434028	<b>DICOM WORKLIST &amp; MPPS</b> Import of patient/examination data from an external RIS (Radiology Information System) /HIS (Hospital Information System) patient management system with DICOM MWL (Modality Worklist) as well as feedback on the examination status with DICOM MPPS (Modality Performed Procedure Step).
1	04453846	<b>CAREPROFILE - C</b> Reduces dose for patients and staff by radiation free collimation. The collimator blades are placed using the last-image-hold for orientation.

Qty	Part No.	Item Description
1	14402070	<b>Fluoroloop - C</b> Storage and review of fluoroscopy sequences for documentation. The maximum storable fluoroscopy time depends on the selected pulse rate.
1	04484114	<b>Security Package</b> Software extension for workplaces enabling enhanced security features including user management and audit trail functionality
1	14436606	<b>VA Kit</b> Second set of documentation for Veterans' Affairs Administration Hospitals in the U.S.
1	14443181	<b>Keyboard, US English</b> PS2 standard keyboard
1	14417233	<b>External DVI interface</b> Standard video output (DVI-D format) at the FLUOROSPOT Compact, for connecting an external recorder to record image information on video recording media.
1	14428635	<b>Trolley with 1 Display, b/w</b> Mobile trolley with one titable flat display in TFT technology (high luminance and extended viewing angle) for live image display in the examination room and radiation indicator.
1	14443215	<b>Ysio Max Ceiling Carriage 3 m</b> Ysio Max Ceiling Carriage 3 m Universal digital radiographic workplace for skeletal radiography of the recumbent, standing or seated patient. High-resolution, permanently installed or wireless detectors as a basis for a fully digital imaging chain with a digital imaging system, an image and control station with application and evaluation programs, and DICOM network connection. Tube assembly support fully motorized in all projection-relevant axes with up to 220 cm transverse travel. OPTITOP 150/40/80 X-ray tube assembly and multileaf collimator with full field and laser line light localizer.
1	14436409	<b>Ceiling rails 4.25m</b> 2 tracks for the ceiling-mounted support with a travel distance up to a maximum of 4.25 meters in longitudinal direction.
1	14407006	<b>Caremax plus HS Integrated</b> CAREMAX plus Dose Area Product (DAP) meter tracks and displays the Dose Area Product (DAP) and/or standardized patient entrance dose.
1	14443192	<b>Bucky Wall Unit for MAX wi-D Left</b> Bucky Wall Unit for MAX wi-D Left Floor-mounted Bucky wall stand with height-adjustable and tiltable detector tray for a MAX wi-D flat detector for digital acquisitions. With IONTOMAT three-field chamber and Bucky frame. Detector tray operated from the left / right side. Vertical height adjustment and detector tilt possible from both sides.
1	14428861	<b>Transparent grid 13/92, Universal</b> Highly selective anti-scatter grid for scattered radiation reduction.
1	14443197	<b>WLAN US</b>
1	14443195	<b>MAX wi-D</b>
1	14436593	<b>MAX wi-D Clip-on Grid 5/85 F115</b> Grid (5/85), f 115 cm Highly selective anti-scatter grid for scattered radiation reduction:

Qty	Part No.	Item Description
		<ul style="list-style-type: none"> <li>- Pb 5/85 (grid ratio 5:1, 85 lines/cm)</li> <li>- Grid focusing for SID 115 cm (45")</li> </ul>
1	14436629	<b>Ortho Stand</b> <p>The Ortho support stabilizes the patient during the Ortho examination, such as long leg and long spine imaging to prevent motion artifacts. It places the patient in an upright position and allows safe movement of the bucky tray during the image acquisition procedure.</p> <p>The package includes:</p> <ul style="list-style-type: none"> <li>- Ortho ruler</li> <li>- Ruler holder</li> <li>- Patient hand grips, left and right</li> <li>- Additional platform for smaller patients (e.g. children) and to make sure that the whole patient body can be mapped - removable</li> </ul>
1	08861002	<b>Patient positioning mattress</b> <p>The radiolucent table pad matches the size of the tabletop and has a heavy-duty soft plastic cover that is easy to clean. The soft cushion allows comfortable patient positioning and repositioning. To prevent the pad from sliding during head-up positions, the straps of patient table pad can be attached to the grip protection rail at the head end.</p>
1	14417389	<b>Customer documentation, English</b> <p>One set of manuals is included with the system. This option is only for customers that require more than one set of system manuals.</p>
1	AXD_INITIAL_24	<b>Initial onsite training 24 hrs</b> <p>Up to 24 hours of on-site clinical education training, scheduled consecutively (Monday - Friday) during standard business hours for a maximum of (4) imaging professionals. Training will cover agenda items on the ASRT approved checklist. Uptime Clinical Education phone support is provided during the warranty period for specified posted hours. This educational offering must be completed (12) months from install end date. If training is not completed within the applicable time period, Siemens obligation to provide the training will expire without refund.</p>
1	AXD_ADD_16	<b>Additional onsite training 16 hours</b> <p>Up to (16) hours of on-site clinical education training, scheduled consecutively (Monday - Friday) during standard business hours for a maximum of (4) imaging professionals. Training will cover agenda items on the ASRT approved checklist if applicable. This educational offering must be completed (12) months from install end date. If training is not completed within the applicable time period, Siemens obligation to provide the training will expire without refund.</p>
1	CID4948	<b>Portable DR Panel Protector(14x17)</b> <p>The unique design of the DR Panel Protector provides an easy way to take weight-bearing x-rays of feet (AP view). The unit is simply placed over the DR panel which is first positioned on the floor. Patients step onto the DR Panel Protector with as much weight as needed to get the desired image. The face plate is made of polycarbonate designed to support patients weighing up to 500 pounds. The face plate is x-ray lucent, allowing the x-rays to pass through the DR Panel Protector with no significant absorption or scattering. The non-slip rubber floor grips keep the DR Panel Protector from slipping on a hard floor. The Panel Protector frame is notched to accommodate the cable connection from the digital DR panel to the host system. One year warranty through Clear Image Devices</p>
1	TIMS061000SP PKGSHC	<b>TIMS 2000 SP Package</b> <p>Includes: TIMS 2000 SP desktop system, 19" monitor, one license of TIMS DICOM Review Software (TDRS), wireless microphone, keyboard, mouse, and associated cables.</p> <p>Description: Provides for the high resolution video and synced audio recording, review, analysis, and archiving to PACS of modified barium swallow studies for speech pathology.</p>

Qty	Part No.	Item Description
		<p>The system is positioned on the desktop of the fluoroscopy room control area. Also includes one license of TIMS DICOM Review Software (TDRS) for the remote review and analysis of studies in the speech pathology office area. This allows for more efficient use of the fluoroscopy room because studies are sent immediately to the TDRS system at the completion of the study. With this workflow, the fluoroscopy room is immediately available for the next patient, allowing more studies per day. Speech pathologists can then edit, review, label, analyze, and archive to PACS studies from the comfort of their own office. Onsite installation and training is included when purchased with a new Siemens' system. When purchased to upgrade an existing system, the customer is responsible for room preparation and cable installation per Foresight Imaging's planning guide. Additional costs may be required for a Siemens service engineer.</p> <p>Warranty: Support &amp; Maintenance for one year. Includes technical support via telephone, email, and online on regular business days from 8:00am to 8:00pm EST, downloadable software updates, and in the event of a hardware failure, a replacement system is sent within two business days.</p>
1	AS10847102	<p><b>Mobile detector holder for Max wi-D</b></p> <p>The versatile holder accommodates computed radiography (CR) cassettes and light portable DR Panels (including the max wi-D detector) with a total weight (including clip-on grid if required) of less than 4.3kg (9.5 lbs).</p> <p>The holder rolls on large locking castors and facilitates examinations in accident and emergency departments, in operating rooms and radiographic rooms. The heavy duty base gives a low center of gravity.</p> <p>Properties:</p> <ul style="list-style-type: none"> <li>• The holder is adjustable for height from floor level to 50 in (measured from its lower edge)</li> <li>• The holder is counterbalanced for easy rising or lowering and can overhang the x-ray or operating table by 24 in.</li> <li>• Supports detectors with a width of 9.6 to 21 in</li> <li>• Maximum detector thickness 1.2 in (including clip-on grid if required)</li> <li>• The holder can be turned &amp; tilted and orientated to suit any examination position</li> <li>• Effective locks keep the holder firmly in place</li> </ul>
1	AXD_RIG_FLU ORO_GOV	<b>Standard Rigging Fluoro</b>
1	AXD_ADDL_RI GGING	<b>Additional Rigging AXD</b>

---

**ALTERNATE PRODUCTS on Quote Nr: 1-LFWL94 Rev. 1**

---

**ALTERNATE PRODUCTS for Luminos Agile Max**

**All items listed below are ALTERNATE PRODUCTS:** *(See Detailed Technical Specifications at end of Proposal.)*

Qty	Part No.	Item Description
1	14428636	<b>DCS 1FT with 1 Display b/w</b> Flexible postioning of two flat displays in TFT technology (high luminance and extended viewing angle) in the examination room to save space on the floor. The display suspension system includes a radiation indicator, is ceiling-mounted, swiveling rotatable and height-adjustable with longitudinal travel  <b>The above item(s) are being quoted as a substitute for the following quoted Part No(s).</b> <ul style="list-style-type: none"><li>• 14428635 - "Trolley with 1 Display, b/w</li></ul>

---

OPTIONS on Quote Nr:

1-LFWL94 Rev. 1

---

## OPTIONS for Luminos Agile Max

All items listed below are **OPTIONS**: *(See Detailed Technical Specifications at end of Proposal.)*

Qty	Part No.	Item Description
1	14443367	<b>SmartOrtho License</b> SmartOrtho is an automated tilting technique for long leg and full spine imaging. Up to 4 single images can be acquired to cover the selected region with the patient in standing or lying position. The images are automatically composed into a single image on the imaging system.
1	14428919	<b>Transparent grid 15/80, F300</b> Highly selective anti-scatter grid for scattered radiation reduction.
1	AXD_ADD_12	<b>Additional onsite training 12 hours</b> Up to (12) hours of on-site clinical education training, scheduled consecutively (Monday - Friday) during standard business hours for a maximum of (4) imaging professionals. Training will cover agenda items on the ASRT approved checklist if applicable. This educational offering must be completed (12) months from install end date. If training is not completed within the applicable time period, Siemens obligation to provide the training will expire without refund.

# Detailed Technical Specifications

## Luminos Agile Max

Part No. / Product	Description
14436607 Luminos Agile Max	<p><b>System configuration</b></p> <p><b>Patient positioning</b> Patient positioning table easily accessible from all sides, tiltable from +90° to max. 20° Trendelenburg position.</p> <ul style="list-style-type: none"> <li>- Table height adjustable from 65 cm to 112 cm (min. fluoro height 82 cm)</li> <li>- Patient positioning tabletop 210 cm x 80 cm with examination range over the entire radiolucent part (193 cm x 53 cm). Motorized longitudinal travel <math>\pm 80</math> cm and transverse travel <math>\pm 17.5</math> cm.</li> <li>- Foot rest attachable at head or foot end, up to a load of 230 kg.</li> <li>- Comprehensive collision protection within the travel range of the image receptor.</li> </ul> <p>Patient weight:</p> <ul style="list-style-type: none"> <li>- Up to 150 kg - No restriction.</li> <li>- 150 kg -182 kg (400 lbs) Patient longitudinal travel <math>\pm 40</math> cm with slow movement, table tilt unrestricted.</li> <li>- Between 181 kg and 230 kg if tabletop is in center position. Table tilt from + 90° to - 15°.</li> <li>- Between 231 kg and 275 kg if tabletop is in center position.</li> </ul> <p>Interventions up to a patient weight of 250 kg (plus 50 kg still possible for CPR).</p> <p><b>Flat detector image receptor</b></p> <ul style="list-style-type: none"> <li>- Distance focus - tabletop 60 cm.</li> <li>- SID from 89 cm to 125 cm.</li> <li>- Primary collimator with rectangular format.</li> <li>- Highly selective anti-scatter grid Pb 15/80, grid ratio 15:1, 80 lines/cm. Grid focusing for SID 100 cm.</li> </ul> <p>Transparent grids improve the image quality by reducing scattered radiation on the film.</p> <p><b>Flat detector:</b> The digital high-resolution dynamic flat detector with integrated removable grid is especially designed to fulfill the requirements of general and interventional applications.</p> <ul style="list-style-type: none"> <li>- Semi-conductor material: Amorphous silicon (a-Si) with CsI-(scintillator)</li> <li>- Size 43 cm x 43 cm</li> <li>- Pixel size: 148 <math>\mu</math>m (6.76 pixel per mm)</li> <li>- Physical size of matrix: 2880 x 2880 (8.3 million pixels)</li> <li>- Size of active matrix: 2840 x 2874 (8.2 million pixels)</li> <li>- Detail resolution: 3.4 LP/mm</li> <li>- Acquisition depth: 16 bits</li> </ul> <p>148 <math>\mu</math>m pixel arrays provide highest spatial resolution and excellent contrast. Fluoroscopy as well as image acquisition are always done in 14-bit gray scale resolution, allowing excellent detail visibility. Acquisition frame rates of up to 8 f/s are possible.</p> <p>Usable input formats:</p> <ul style="list-style-type: none"> <li>- Overview: 43 cm x 43 cm; diagonal 60 cm.</li> <li>- Zoom 1: 30 cm x 30 cm; diagonal 42 cm.</li> </ul>

Part No. / Product	Description
<p><b>(Continued)</b>  <b>14436607</b>  <b>Luminos Agile Max</b></p>	<ul style="list-style-type: none"> <li>- Zoom 2: 22 cm x 22 cm; diagonal 32 cm.</li> <li>- Zoom 3: 15.5 cm x 15.5 cm; diagonal 21 cm.</li> </ul> <p><b>Operation</b>  Ergonomic tableside user interface on the system, comfortable for both left- and right-handed users. Integrated system operation and total tableside examination control.</p> <ul style="list-style-type: none"> <li>- OPTI Grip for fast and easy movement of the flat detector image receptor and single-handed system operation.</li> <li>- Touch user interface for interactive control of all system functions.</li> </ul> <p>Footswitch FL/R:  The footswitch combination with two separate switch pedals provides ergonomic work support. During fluoroscopy-guided patient examinations, the user's hands remain free for the examination.  Fluoroscopy is switched on and off by foot. During fluoroscopy, targeted radiographic acquisitions can be released by foot, as well.</p> <p><b>Accessories included in basic version</b></p> <ul style="list-style-type: none"> <li>- Hand grip, front, axial adjustment.</li> <li>- Hand grip rail, back, axial adjustment.</li> <li>- Grip protection rail, head end, removable.</li> <li>- Shoulder supports (1 pair), three-dimensional adjustment.</li> <li>- Foot board, axial adjustment by 44 cm, can be used at head end or foot end.</li> <li>- Protective film for fluids</li> <li>- Radiation protection, complete</li> <li>- Counterweight for radiation protection</li> </ul> <p><b>OPTITOP 150/40/80HC-100 undertable X-ray tube assembly</b>  X-ray tube OPTITOP 150/40/80 HC-100:</p> <ul style="list-style-type: none"> <li>- Single-track, dual-focus rotating anode tube with compound anode (rhenium-tungsten, molybdenum, graphite)</li> <li>- High heat storage capacity and high thermal load capacity for small focal spots. Integrated overpressure safety device in the tube protective housing.</li> <li>- 150 kV nominal voltage acc. to IEC 60613.</li> <li>- Nominal radiographic anode input power acc. to IEC 60613 (focal spot nominal values acc. to IEC 60336):  47 kW: small focus 0.6  85 kW: large focus 1.0</li> <li>- Anode speed <math>\geq</math> 8,500 rpm, anode angle 12°.</li> <li>- Heat storage capacity of the anode 580 kJ (783 kHU) acc. to IEC 60613.</li> <li>- Nominal radiographic anode input power acc. to IEC 60613 (focal spot nominal values acc. to IEC 60336)</li> </ul> <p><b>Generator Polydoros</b>  High-frequency X-ray generator with multipulse voltage waveform for diagnostic acquisition techniques. The multipulse voltage waveform enables high data accuracy, precise reproducibility and short exposure times.</p> <ul style="list-style-type: none"> <li>- Rating:  65 kW at 100 kV acc. to IEC 60601-2-7.  max. 800 mA at 79 kV</li> <li>- Tube voltage 40 kV to 150 kV.</li> <li>- Fluoroscopy, 450 W from 40 kV/0.2 mA to 110 kV/23 mA, pulsed fluoroscopy</li> <li>- Shortest exposure time 1 ms (with IONTOMAT automatic exposure control). Free selection of the exposure data in 3- and 2-point technique.</li> <li>- 1-point technique with continuously falling load (with IONTOMAT).</li> <li>- Plani-IONTOMAT for tomography (option).</li> <li>- CAREMATIC system for 0-point technique.</li> </ul>



Part No. / Product	Description
<p><b>(Continued)</b>  <b>14436607</b>  <b>Luminos Agile Max</b></p>	<ul style="list-style-type: none"> <li>- Time and mAs post indication (with IONTOMAT operation).</li> <li>- Fluoro mean value indication.</li> <li>- Organ programs can be edited by the user.</li> <li>- Tube load computer with opto-acoustic warning indicator.</li> </ul> <p>Power connection:  3 phase current: 400 V (-15, +10%); 50/60 Hz.  Option: 440/480 V.</p> <p><b>Multileaf collimator</b>  Undertable multileaf collimator with rectangular collimation for automatic format collimation. Motor-driven Cu prefilters.</p> <p><b>CAREMAX</b>  Electronic unit with KermaX-Plus, a measurement chamber integrated into the collimator housing for acquisition and fluoro systems to record the dose area product and/or standardized patient entry dose.</p> <p><b>19" Flat display</b></p> <ul style="list-style-type: none"> <li>- TFT flat-screen display for live image display in the control room.</li> <li>- Flicker-free and distortion-free image display.</li> <li>- Screen size: 19"(48 cm).</li> <li>- Resolution: 1280 x 1024 pixels.</li> <li>- Maximum brightness: <math>\leq 800 \text{ cd/m}^2</math>.</li> </ul> <p><b>Imaging system</b>  High-resolution digital imaging system with innovative image display, DICOM network connection and syngo-like user interface. It was optimized for general fluoro and multi-functional workplaces.</p> <p>With the FLUOROSPOT Compact both single acquisitions and series from 0.5 to 8 images per second in 1440<sup>2</sup> matrix can be made and reviewed dynamically.  The digital acquisition process produces continuously high-quality images for all fluoroscopy-guided contrast medium examinations, skeletal acquisitions and interventional procedures.</p> <p><u>Operating modes:</u></p> <ul style="list-style-type: none"> <li>- Digital radiography (DR) with up to 2, 880<sup>2</sup>/12-bit matrix.</li> <li>- CAREVISION: Pulsed fluoroscopy with selectable pulse frequencies 3 f/s, 7.5 f/s, 10 f/s, or 15 f/s in 1024<sup>2</sup>/12-bit matrix.</li> <li>- High-speed fluoro at 30 f/s in Zoom 2 (22 cm x 22 cm) and Zoom 0 (43 cm x 43 cm)  Pulsed fluoroscopy is especially suitable for time-intensive examinations to reduce the radiation dose for physician, staff, and patient.</li> <li>- Display and storage of the last fluoroscopic image after switching off radiation (Last Image Hold).</li> <li>- Single image and serial mode up to 1,440<sup>2</sup>/12-bit or 1,024<sup>2</sup>/12-bit matrix (depending on zoom).</li> <li>- Serial mode (max. 8 f/s) with variable frame rate in three steps (max. 1,440<sup>2</sup>/12-bit matrix).</li> <li>- DDO (Dynamic Density Optimization) for on-line harmonization of native series and single images.</li> <li>- Live auto-windowing during fluoroscopy</li> <li>- Live auto-shutter during fluoroscopy</li> </ul> <p><u>Image processing:</u></p> <ul style="list-style-type: none"> <li>- Real-time edge enhancement, positive/negative image display, windowing, contrast/brightness, electronic display (shutter), image shift (roaming), vertical and horizontal image inversion, and magnifying glass and zoom functions.</li> <li>- DiamondView Plus: multi-scaling procedure for image post-processing with high detail contrast and reduced noise. Filter size and strength are weighted differently and are used for adaptation to the overall image content. DiamondView Plus enhances the signal exploitation of the dynamic range and improves the organ-</li> </ul>

Part No. / Product	Description
<b>(Continued)</b> <b>14436607</b> <b>Luminos Agile Max</b>	<p>specific detail contrast (soft tissue and bone). Can be preset in the organ program, and can subsequently be selected or deselected.</p> <ul style="list-style-type: none"> <li>- Storing of single images as reference images also during fluoroscopy.</li> <li>- Quantification: angle/length measurement, automatic and/or manual calibration.</li> <li>- Text functions: User-defined image annotation, free annotation or using text components, comment line for the image, R/L display, image report and image quality graphics.</li> <li>- Image gallery for harmonization (visualization support for image post-processing)</li> </ul> <p><u>Image display:</u> Image review and display in 100 Hz progressive display (1,024<sup>2</sup> matrix) through high-resolution, flicker-free flat-screen displays.</p> <ul style="list-style-type: none"> <li>- Screen layout with 4, 9, 16, 25, or 36 images of an examination.</li> <li>- Display of R/L marks.</li> </ul> <p><u>Image storage capacity:</u> 50,000 images for permanent storage in 1k/12-bit matrix and 2,000 images for permanent storage in 2840 x 2880 matrix.</p> <p><b>DVD / CD burner (DICOM)</b> DVD drive for automatic digital image storage on CD/DVD for offline data exchange in DICOM, TIFF, and AVI format.</p> <p><b>Connectivity</b></p> <ul style="list-style-type: none"> <li>- DICOM Send: Digital, unidirectional image transfer of single images or complete folders to a network in DICOM format.</li> <li>- DICOM Print: Provision of DICOM Print service for connection to a laser camera or a network printer (postscript-capable).</li> <li>- DICOM Storage Commitment (StC): The network/archive sends a receipt acknowledgment for images/folders to the image system in DICOM format.</li> <li>- DICOM Dose structured report: Enables transfer of dose report to archiving or viewing station.</li> </ul> <p><b>Note concerning DICOM interface(s)</b> For diagnostic purposes, only hardcopy cameras/laser printers explicitly approved for this system may be used.</p> <p>The description in the "DICOM Conformance Statement" downloadable from the Internet is exclusively binding for the functionality of the DICOM interface(s).</p> <p>Functionalities across system borders with/between partner systems require explicit validation, since the interpretation of the interface by the partner/target system is not part of the product's responsibility.</p> <p>A modification of the interface that might be required is not included in the offer; e.g. for the rare case, that available configurations are not sufficient. With regard to expenses for interface configurations that might be required, the agreements on maintenance/service of the product apply.</p> <p><b>syngo Remote Assist</b> <i>syngo</i> Remote Assist is a standalone service option. With <i>syngo</i> Remote Assist, Siemens uses a secure broadband VPN connection (VPN = virtual private network) to establish a connection to your Siemens imaging console in order to offer you direct, real-time support and training. This seamless and simultaneous virtual interaction will contribute to improvements in image quality and optimization of system use.</p> <p><b>Siemens Remote Service</b> System Management software package to support Siemens Remote Service (SRS) with the following functions:</p> <ul style="list-style-type: none"> <li>- Basic package Siemens Remote Service for Diagnostics and Repair, Quality Assurance and Software Maintenance.</li> <li>- System remote configuration, e.g. adding of a DICOM node.</li> </ul>

Part No. / Product	Description
<b>(Continued)</b> <b>14436607</b> <b>Luminos Agile Max</b>	<ul style="list-style-type: none"> <li>- Early warning system to secure system operation.</li> <li>- The functions are made available in accordance with the maintenance contract package.</li> </ul> <p>Prerequisite for the early warning system is a permanent connection to the system via LAN and router. It is the project manager's task to make this available on-site.</p> <p><b>Power connection for the entire X-ray system</b>  Basic version: 3/N/PE ~400 V (<math>\pm 10\%</math>) for 50/60 Hz line frequency.  Option: 3/N/PE ~440/480 V (<math>\pm 10\%</math>) at 50/60 Hz via line adaptation transformer.</p> <p><b>Customer Care. Life - the customer care solution by Siemens Healthcare</b>  From the moment you purchase your Siemens system you will benefit from many services that are offered by Customer Care. Life* offers, e.g.:</p> <ul style="list-style-type: none"> <li>- initial application training,</li> <li>- interactive e-learning for various applications,</li> <li>- free customer magazines,</li> <li>- arrangements for clinical training via a global network,</li> <li>- and free trial licenses</li> </ul> <p>You will find detailed information on our e-learning program and further details on general Customer Care. Life services on the internet.</p> <p>* Not all services of the Customer Care. Life offerings are necessarily available for all systems.</p> <p>Components for basic configuration are described in the following.</p>
<b>14417363</b> <b>Upgrade to 80 kW</b> <b>PLF</b>	Performance data : <ul style="list-style-type: none"> <li>- 80 kW at 100 kV (IEC 60601-2-7), 1,000 mA at 80 kV.</li> <li>- Fluoroscopy, 450 W from 40 kV/0.2 mA to 110 kV/23 mA, pulsed fluoroscopy</li> </ul>
<b>14428860</b> <b>Transparent grid</b> <b>13/92, F115</b>	Technical details: <ul style="list-style-type: none"> <li>- Grid ratio 13:1, 92 lines/cm</li> <li>- Grid focusing for source-image-distance (SID) of 115 cm (45")</li> </ul>
<b>04434028</b> <b>DICOM WORKLIST &amp; MPPS</b>	<p><b>DICOM MWL (Modality Worklist):</b>  Import of patient/examination data from an external RIS/HIS patient management system.</p> <p><b>DICOM MPPS (Modality Performed Procedure Step):</b>  Sending of dose data, patient data, and examination data to an external RIS/HIS patient management system.</p> <p><b>Note concerning DICOM interface(s)</b>  The description in the "DICOM Conformance Statement" downloadable from the Internet is exclusively binding for the functionality of the DICOM interface(s).</p> <p>Functionalities across system borders with/between partner systems require explicit validation, since the interpretation of the interface by the partner/target system is not part of the product's responsibility.</p> <p>A modification of the interface that might be required is not included in the offer; e.g. for the rare case, that available configurations are not sufficient.  With regard to expenses for interface configurations that might be required, the agreements on maintenance/service of the product apply.</p>
<b>04484114</b> <b>Security Package</b>	<ul style="list-style-type: none"> <li>- This SW license enables the Workplace System to support enhanced user and system management,</li> </ul>

Part No. / Product	Description
<b>(Continued)</b> <b>04484114</b> <b>Security Package</b>	<p>including: user authentication to prohibit unauthorized access.</p> <ul style="list-style-type: none"> <li>- privileges to define user/role based functionality.</li> <li>- permissions to control data access.</li> <li>- audit trails to log system and data access.</li> </ul>
<b>14428635</b> <b>Trolley with 1</b> <b>Display, b/w</b>	<p>Technical details of display:</p> <ul style="list-style-type: none"> <li>- Screen size 19" (48 cm)</li> <li>- Resolution: 1280 x 1024 pixels</li> <li>- Maximum brightness (typical): <math>\leq 800 \text{ cd/m}^2</math></li> <li>- Flicker-free and distortion-free image display</li> <li>-</li> </ul>
<b>14443215</b> <b>Ysio Max Ceiling</b> <b>Carriage 3 m</b>	<p><b>System Configuration</b></p> <p>Ysio Max is a universal digital radiographic workplace with various flat detectors (MAX wi-D, MAX static) for image acquisition.</p> <p>The Ysio Max digital workplace is especially suited for a high patient throughput. As a universal workplace, the system is primarily used in X-ray departments of hospitals, in radiological and partly radiological offices with high patient throughput and standardized acquisition technology.</p> <p>Basic system components:</p> <ul style="list-style-type: none"> <li>- A ceiling-mounted tube assembly support with X-ray tube assembly and motorized multileaf collimator.</li> <li>- An imaging and control station with application and evaluation programs, as well as DICOM system interfaces.</li> <li>- CD/DVD drive for digital image storage on CD-R/DVD for offline data exchange in DICOM format.</li> </ul> <p><b>Tube assembly support</b> with X-ray tube assembly and motorized collimator.</p> <p>All projection-relevant tube assembly positions can be manually adjusted with handles symmetrically mounted to the tube assembly collimator unit.</p> <p>The ceiling-mounted tube assembly support can be adjusted in 3 axes for longitudinal, transverse, and height adjustment (x, y, and z-axes).</p> <ul style="list-style-type: none"> <li>- Horizontal travel range in longitudinal direction 346 cm.</li> <li>- Horizontal travel range in transverse direction 220 cm.</li> <li>- Vertical lift 180 cm.</li> </ul> <p>In 2 further axes (<math>\alpha</math>- and <math>\beta</math>-axes) the tube assembly collimator unit can be manually adjusted for oblique acquisitions of the recumbent patient, or for horizontal, oblique, or lateral acquisitions on the portable detector, or for free bedside acquisitions.</p> <ul style="list-style-type: none"> <li>- Rotation around the vertical axis of the ceiling-mounted support from <math>+154^\circ</math> to <math>-182^\circ</math>. Lock-in positions every <math>90^\circ</math>.</li> <li>- Rotation around the horizontal axis of the tube assembly support arm <math>\pm 140^\circ</math>. Lock-in positions at <math>0^\circ</math> and <math>\pm 90^\circ</math>.</li> </ul> <p>X-ray tube assembly OPTITOP 150/40/80 HC-100: Single-track dual-focus rotating anode tube with compound anode (rhenium-tungsten, molybdenum, graphite), with high heat storage capacity and high load capacity for small focal spots. Integrated overpressure safety device in the tube protective housing.</p> <ul style="list-style-type: none"> <li>- 150 kV nominal voltage acc. to IEC 613.</li> <li>- Nominal power (focal spot nominal values acc. to IEC 336): 40 kW: small focus 0.6 80 kW: large focus 1.0</li> </ul>

Part No. / Product	Description
<p><b>(Continued)</b>  <b>14443215</b>  <b>Ysio Max Ceiling</b>  <b>Carriage 3 m</b></p>	<ul style="list-style-type: none"> <li>- Anode speed <math>\geq 8,500</math> RPM, anode angle <math>12^\circ</math>.</li> <li>- Heat storage capacity of the anode 580 kJ (783 kWh) acc. to IEC 613.</li> <li>- Total filtration (IEC 601-1-3) <math>\geq 2.5</math> mm Al equiv.</li> </ul> <p>Multileaf collimator:  With full field and laser line light localizer. Rectangular collimation, manual and motorized, via organ programs.</p> <ul style="list-style-type: none"> <li>- Multileaf collimator rotatable by <math>\pm 45^\circ</math> around the center beam axis, e.g. for correct positioning of objects.</li> <li>- A tape measure is integrated to check the focus-to-object distance.</li> <li>- To improve radiation quality through dose reduction of the soft radiation parts, Cu filters (0.1 Cu; 0.2 Cu and 0.3 Cu) are inserted into the primary beam projection, depending on the organ program selected. They can also be selected manually.</li> </ul> <p>Option:  A measuring chamber for the dose area product can be integrated into the multileaf collimator.</p> <p><b>Controls and displays</b>  The control elements at the tube assembly and the multileaf collimator are ergonomically arranged for single-handed operation.</p> <p>Controls and displays at the tube assembly support (MAXTouch):  Multifunctional control display with color touchscreen for adaptation of acquisition parameters directly in the examination room.</p> <p>Displays include:</p> <ul style="list-style-type: none"> <li>- The collimation size of the acquisition field (in cm x cm).</li> <li>- The selected SID.</li> <li>- The selected Cu additional filters.</li> <li>- Rotation from the 0-position.</li> <li>- Tube assembly and detector centering.</li> <li>- Operating states such as "ACSS/Manual", "Ready", "Selected", etc.</li> <li>- Current detector angle (MaxAlign feature) – to eliminate the need to guess the tube angle and to protect the patient by reducing repeat exposures. Available with MAX wi-D and MAX mini.</li> </ul> <p>The display follows the tube assembly orientation.</p> <p>The following functions can be set manually at the multileaf collimator:</p> <ul style="list-style-type: none"> <li>- Full field light localizer with timer for optical display of the collimated acquisition format and an optionally coverable laser line light localizer.</li> <li>- The collimation of the acquisition format set last can be retrieved via a memory button.</li> <li>- The rectangular collimation of the radiation field is pre-defined through the organ program and can be set manually by means of two dials.</li> <li>- The motorized insertion of the Cu additional filters is controlled via the organ program, but can also be selected freely.</li> </ul> <p><b>Imaging and control station (syngo FLC)</b>  The entire control and communication of the radiography system incl. digital image processing takes place from a central operating site - the imaging and control station.</p> <p>It includes:</p> <ul style="list-style-type: none"> <li>- A high-end PC imaging system, based on Windows 7 with syngo user interface.  Storage of original data 14 bit.  Storage of image data 12 bit.  Storage capacity approx. 10,000 images.</li> <li>- Keyboard and mouse.</li> <li>- One 19" color flat-screen as control display or diagnostic display.</li> </ul>

Part No. / Product	Description
<p><b>(Continued)</b>  <b>14443215</b>  <b>Ysio Max Ceiling</b>  <b>Carriage 3 m</b></p>	<ul style="list-style-type: none"> <li>- Manual button for exposure release.</li> </ul> <p><b>Functions of the imaging and control station</b></p> <p>Patient and study administration:</p> <ul style="list-style-type: none"> <li>- Importing of patient lists and examinations from the HIS/RIS</li> <li>- Manual patient registration</li> <li>- Patient, study, and image data management</li> <li>- Configuration functions</li> </ul> <p>Acquisition and postprocessing:</p> <ul style="list-style-type: none"> <li>- Organ program selection and configuration</li> <li>- Selection of generator and diaphragm parameters.</li> <li>- Parameterization of image preprocessing: enhancement, harmonization, edge enhancement, and look-up tables (LUT)</li> <li>- Display of current acquisition between 1.5 and 3.5 seconds (preview); complete image in 3.5 to 6 seconds max. depending on detector type</li> <li>- Display of image markers (L/R, a.p./p.a.)</li> <li>- DiamondView Plus: multi-scaling procedure for image post-processing with high detail contrast and reduced noise</li> </ul> <p>DiamondView is a multi-scale procedure, i.e. filter size and strength are weighted differently and are used for adaptation to the overall image content.</p> <ul style="list-style-type: none"> <li>- DiamondView enhances the signal exploitation of the dynamic range and improves the organ-specific detail contrast (soft tissue and bone).</li> <li>- DiamondView can be selected via the "Pre-processing card".</li> <li>- By entering "0", the image can be displayed without DiamondView.</li> </ul> <p>Image processing functions:</p> <ul style="list-style-type: none"> <li>- Image rotation</li> <li>- Horizontal/vertical image mirroring</li> <li>- Image zoom</li> <li>- Pan</li> <li>- Windowing</li> <li>- Filters for edge enhancement and noise reduction</li> </ul> <p>Image documentation and archiving:</p> <ul style="list-style-type: none"> <li>- Image transfer to the network</li> <li>- Automatic, user-configurable data distribution (DICOM Send, see also system interfaces DICOM)</li> <li>- Automatic filming with virtual film sheet (DICOM Print, see also system interfaces DICOM)</li> <li>- Image data export (12 bit) on CD/DVD</li> </ul> <p><b>Workflow</b></p> <p>Routine workflows are largely automated.</p> <ul style="list-style-type: none"> <li>- Prior to exposure the patient data is transferred via the patient management system (HIS/RIS: option) or entered through the control console. The exposure parameters are selected through the organ programs.</li> <li>- Then the patient or the acquisition system is positioned and exposure is released.</li> <li>- The exposure released at the central system control is read out within a few seconds by the detector. It is displayed at the control display for orientation and made available in DICOM format at the imaging system output for sending e.g. to reporting workstations, image networks, laser cameras, etc.</li> <li>- Clinical Assurance Program (CAP): Collection of deleted images, studies and patient data, including evaluation capabilities.</li> </ul>

Part No. / Product	Description
<p><b>(Continued)</b>  <b>14443215</b>  <b>Ysio Max Ceiling</b>  <b>Carriage 3 m</b></p>	<p>Password protection:  System access protected by password.</p> <p>Option:  Security Package: SW option with enhanced security features such as User Management and Audit Trail function (if offered, see text of the corresponding components).</p> <p><b>DICOM system interfaces</b></p> <ul style="list-style-type: none"> <li>- DICOM Send: Sending of images into the DICOM network.  The DICOM Send function enables fully automatic transfer of generated image data to a DICOM archive or a DICOM workstation. The user can perform his examinations without interruption while the system fully automatically transfers the images to the archive. This image data transfer takes place entirely in the background and thus does not affect acquisitions performed at the same time.</li> <li>- DICOM Storage Commitment (StC): Feedback from the image archive.  The DICOM StC function automatically gives feedback on whether the generated image data were successfully transferred. This way the user can be sure that the acquisitions stored locally in the imaging system can be deleted.</li> <li>- DICOM Print: Printing of images by means of a virtual filmsheet on a DICOM laser camera.  Selecting "Auto-Print" automatically forwards the images stored in the virtual filmsheet to the laser camera. This optimizes the workflow, eliminating the need for user interaction. In addition, a specific layout can be configured on the virtual filmsheet, which the user can review and edit on the monitor at any time. As a result, printing is only required after the layout has been optimized on the monitor, saving time and costs.</li> </ul> <p>Options:</p> <ul style="list-style-type: none"> <li>- DICOM Modality Worklist/MPPS (if offered, see tender further down).</li> <li>- DICOM Query/Retrieve (if offered, see tender further down).</li> </ul> <p><b>Note concerning DICOM interface(s)</b>  For diagnostic purposes, only hardcopy cameras/laser printers explicitly approved for this system may be used.</p> <p>The description in the "DICOM Conformance Statement" downloadable from the Internet is exclusively binding for the functionality of the DICOM interface(s).</p> <p>Functionalities across interfaces with/between partner systems require explicit validation, since the interpretation of the interface by the partner/target system is not part of the product's responsibility.</p> <p>A modification of the interface that might be required is not included in the offer; e.g. for the rare case that available configurations are not sufficient. With regard to expenses for interface configurations that might be required, the agreements on maintenance/service of the product apply.</p> <p><b>syngo Remote Assist</b>  syngo Remote Assist is a standalone service option.  With syngo Remote Assist, Siemens uses a secure broadband VPN connection (VPN = virtual private network) to establish a connection to your Siemens imaging console in order to offer you direct, real-time support and training. This seamless and simultaneous virtual interaction will contribute to improvements in image quality and optimization of system use.</p> <p><b>Siemens Remote Service</b>  Prepared for optional Siemens Remote Service SRS (during warranty period, subsequently with service contract):</p> <ul style="list-style-type: none"> <li>- Hardware and software remote diagnosis.</li> <li>- System remote configuration, e.g. adding of a DICOM node.</li> <li>- Early warning system to secure system operation.</li> <li>- Functions according to the selected maintenance package.</li> </ul> <p><b>Customer Care. Life - the customer care solution by Siemens Healthcare</b>  From the moment you purchase your Siemens system you will benefit from many services that are offered by</p>

Part No. / Product	Description
<p><b>(Continued)</b>  <b>14443215</b>  <b>Ysio Max Ceiling</b>  <b>Carriage 3 m</b></p>	<p>"Customer Care. Life"*, e.g.:</p> <ul style="list-style-type: none"> <li>- initial application training,</li> <li>- interactive e-learning for various applications,</li> <li>- free customer magazines,</li> <li>- arrangements for clinical training via a global network,</li> <li>- and free trial licenses</li> </ul> <p>You will find detailed information on our e-learning program and further details on general "Customer Care. Life" services on the Internet.</p> <p>* "Customer Care. Life" offerings are not necessarily available to the full extent for all systems.</p>
<p><b>14407006</b>  <b>Caremax plus HS</b>  <b>Integrated</b></p>	<p>CAREMAX plus Dose Area Product (DAP) meter is connected to the collimator via CAREMAX adapter cable. The Dose Area Product (DAP) is being displayed on the system control console TUI and recorded in the exam protocol.</p>
<p><b>14443192</b>  <b>Bucky Wall Unit for</b>  <b>MAX wi-D Left</b></p>	<p><b>System Configuration</b>  The Bucky wall stand is a floor-mounted, stand-alone or wall-mountable grid acquisition system with a height-adjustable and tiltable detector tray with tray support and the ability to insert a MAX wi-D flat detector as the digital image acquisition system. The pullout for the opening is left/right.</p> <p>It is especially suited for acquisitions of skeletal radiography of the standing and seated patient:</p> <ul style="list-style-type: none"> <li>- Orthopedic diagnostics.</li> <li>- Thorax and general diagnostics.</li> <li>- Trauma and ER diagnostics.</li> </ul> <p>With this Bucky wall stand, more profound diagnostic requirements for acquisitions of thorax (lungs), abdomen, pelvis, spine, skull and extremities are met.</p> <p>The basic configuration consists of a radiography system with a vertically positioned and tiltable detector Bucky for horizontal, oblique or lateral patient acquisitions.  The additional tilting range of the detector Bucky extends the diagnostically relevant acquisition projections.</p> <ul style="list-style-type: none"> <li>- Vertical height adjustment of the counter-balanced, easily movable detector Bucky from detector center approx. 27 cm to 172 cm above floor: Operation possible from both sides.</li> <li>- Tilting range between 0° and +90°, and up to -20° continuously around the horizontal axis; lock-in position at 0°. Operation possible from both sides.</li> </ul> <p><b>Detector Bucky</b>  The detector Bucky with single-handed operation includes a IONTOMAT three-field chamber for automatic exposure control (incl. three-field templates) and a device for symmetric positioning of the flat detector.</p> <ul style="list-style-type: none"> <li>- Front plate – detector distance ≤42 mm.</li> <li>- Radiation absorption of the front plate ≤0.5 mm Al.</li> <li>- A stationary, exchangeable transparent grid for scattered radiation reduction; Pb 13/92. Optionally for <u>SID</u> 115 cm and/or 180 cm, or Universal Grid with a field from 115 to 180 cm (see tender further down).</li> </ul> <p><b>Detector tray to hold a MAX wi-D detector</b></p> <ul style="list-style-type: none"> <li>- The MAX wi-D detector is supplied with power in the tray and the detector's battery is charged.</li> </ul> <p>The MAX wi-D detector can be inserted in portrait or landscape position. In landscape, the detector is set up for lung acquisitions at the upper edge.  The MAX wi-D detector is detected in the tray.  CR cassettes can also be used directly in the tray.</p> <p><b>Accessories</b></p>



Part No. / Product	Description
<b>(Continued)</b> <b>14443192</b> <b>Bucky Wall Unit for</b> <b>MAX wi-D Left</b>	Scope of delivery: <ul style="list-style-type: none"> <li>- Lateral patient handles for optimum patient positioning, e.g. during PA thorax exposures.</li> <li>- Patient overhead handle, swiveling around the horizontal axis, for optimal patient positioning for lateral acquisitions.</li> </ul>
<b>14428861</b> <b>Transparent grid</b> <b>13/92, Universal</b>	Technical details: <ul style="list-style-type: none"> <li>- Grid ratio 13:1, 92 lines/cm</li> <li>- Grid focusing for source-image distance (SID) of 140 cm (55")</li> <li>- Working range (SID) 115 cm to 180 cm (45" to 71")</li> </ul>
<b>14443195</b> <b>MAX wi-D</b>	Technical details: <ul style="list-style-type: none"> <li>- CsI scintillator, amorphous silicon (a-Si).</li> <li>- Detector acquisition matrix approx. 2872 x 2354</li> <li>- Pixel size 148 µm</li> <li>- Acquisition depth (gray scales) 16 bit</li> <li>- Acquisition formats up to 34,9 cm x 42,5 cm (13,7" x 17")</li> <li>- Thickness 19mm</li> <li>- Data transfer via W-LAN</li> <li>- Operation time: <ul style="list-style-type: none"> <li>- min. 3.5 hours during regular utilization</li> <li>- min. 6 hours in standby mode</li> </ul> </li> <li>- Detector weight 3 kg</li> <li>- Max. load 150 kg (patient lying down) and 100 kg (patient standing)</li> </ul>
<b>14436593</b> <b>MAX wi-D Clip-on</b> <b>Grid 5/85 F115</b>	Technical details: <ul style="list-style-type: none"> <li>- Grid ratio 5:1, 85 lines/cm</li> <li>- Grid focusing for source-image-distance (SID) of 115 cm (45")</li> <li>- Dimensions (W x H x D): 472.1 mm x 410.1 mm x 28.4 mm (18.58"x 16.14"x 1.1")</li> <li>- Weight: 1.1 kg (2.4 lbs)</li> </ul>
<b>14436629</b> <b>Ortho Stand</b>	Technical details: <ul style="list-style-type: none"> <li>- Dimensions (depth x width x height): 75 cm x 75 cm x 202 cm (30" x 30" x 80")</li> <li>- Weight: 85 kg ( 187 lbs)</li> <li>- Maximum patient weight capacity: 180 kg (397 lbs)</li> <li>- Patient height: Up to 190 cm / 74" (standing)</li> </ul>
<b>08861002</b> <b>Patient positioning</b> <b>mattress</b>	Technical details: <ul style="list-style-type: none"> <li>- Length: 198 cm (78")</li> <li>- Width: 66 cm (26") (of which 53.5 cm is padded)</li> <li>- Thickness: 2.5 cm (1")</li> <li>- Weight: 2.7 kg (5,9")</li> </ul>
<b>14428636</b> <b>DCS 1FT with 1</b> <b>Display b/w</b> <b>(Alternate)</b>	Technical details: <ul style="list-style-type: none"> <li>- Screen size 19" (48 cm)</li> <li>- Resolution: 1280 x 1024 pixels</li> <li>- Maximum brightness (typ.): &lt;=800 cd/m<sup>2</sup></li> <li>- Flicker-free and distortion-free image display</li> </ul>

Part No. / Product	Description
<b>14443367</b> <b>SmartOrtho License</b> <b>(Optional)</b>	<p><b>SW license for Ortho function</b></p> <ul style="list-style-type: none"> <li>- For systems with bucky wall stand: Ability to acquire up to 4 images of the legs or spine in sequence on the Bucky wall stand using a MAX wi-D, MAX static detector</li> </ul> <p><u>Only for Ysio:</u> Ability to acquire up to 3 images at the patient table.</p> <p><u>Only for Multitom Rax:</u> Ability to acquire up to 4 images of the legs or spine with RAX detector (patient in standing position or lying on the table)</p> <p><b>Spine Composing</b> Spine Composing takes individually acquired digital radiographic images of the spine and composes them into an overall image.</p> <p><u>The main functions are:</u></p> <ul style="list-style-type: none"> <li>- automatic composing of digital radiographs into an overall image</li> <li>- standard image post-processing functions are available</li> </ul> <p><b>Ortho-Leg Composing</b> Ortho-Leg Composing takes individually acquired digital radiographic images of the legs and</p> <p><u>The main functions are:</u></p> <ul style="list-style-type: none"> <li>- automatic composing of digital radiographs into an overall image</li> <li>- standard image post-processing functions are available composes them into an overall image.</li> </ul>
<b>14428919</b> <b>Transparent grid</b> <b>15/80, F300</b> <b>(Optional)</b>	<p>Technical details:</p> <ul style="list-style-type: none"> <li>- Grid ratio 15:1, 80 lines/cm</li> <li>- Siemens lead/fibre technology</li> <li>- Grid focussing for SID (source image distance) of 300 cm (118")</li> </ul>

# LUMINOS AGILE / AGILE MAX SYSTEM

## TYPICAL ROOM PLAN

---




The intended use for this Cut Sheet is to communicate the spatial requirements as well as the basic architectural, electrical, structural, and mechanical requirements for this piece of imaging equipment. The information provided in this document is for reference only, during the pre-planning stage, and therefore does not contain any site specific detailed requirements. This information is subject to change without notice. Federal, state and/or local requirements may impact the final placement of the components. It is the customer's responsibility to ensure that the final layout and placement of the equipment complies with all applicable requirements.

FOR REFERENCE ONLY,  
NOT FOR CONSTRUCTION.

# LUMINOS AGILE / AGILE MAX SYSTEM

## TYPICAL ROOM PLAN

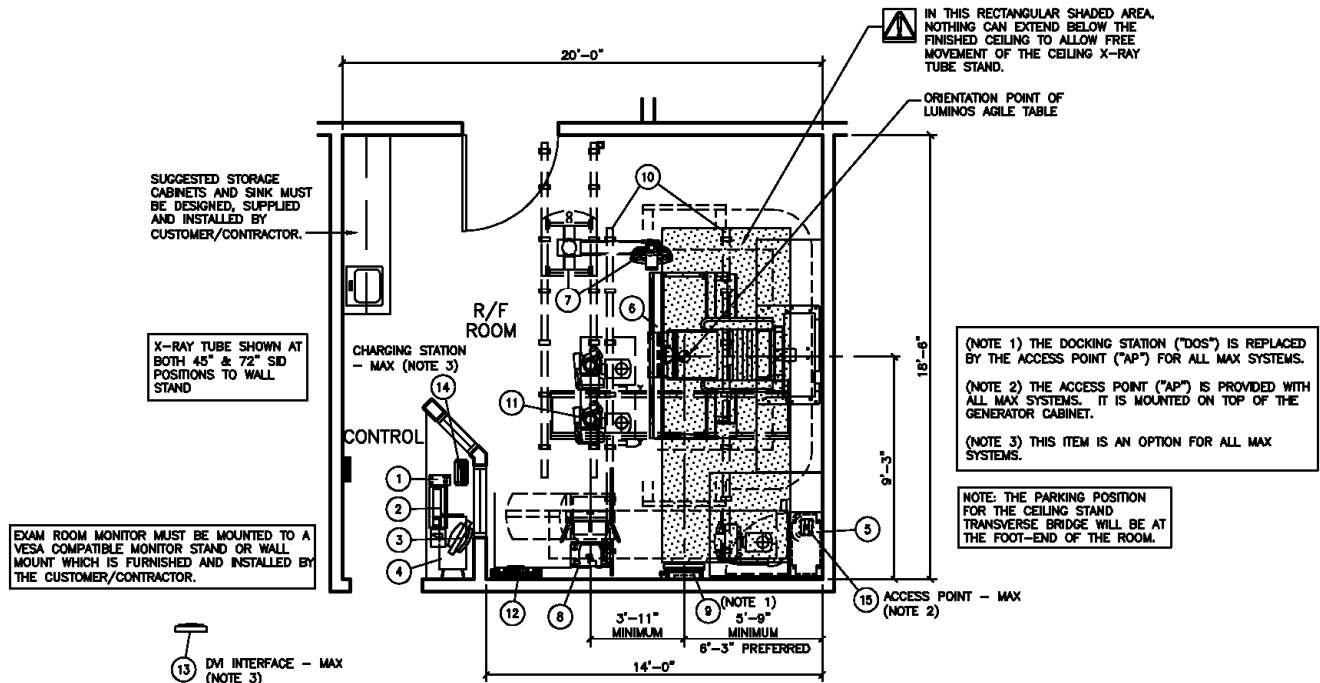

## TYPICAL PLAN

SCALE: 1/8" = 1'-0"

EQUIPMENT LEGEND								
NO	DESCRIPTION	SMS SYM	WEIGHT (LBS)	BTU/HR TO AIR	DIMENSIONS (INCHES)			REMARKS
					W	D	H	
①	CONTROL ROOM MODULE	Ⓢ	2	—	4 3/4	10 1/16	3 3/16	ON CUSTOMER'S COUNTER
②	IMAGING SYSTEM - KEYBOARD AND MOUSE	Ⓢ	—	—	—	—	—	ON CUSTOMER'S COUNTER
③	B/W FLAT SCREEN CONTROL ROOM MONITOR	Ⓢ	—	—	—	—	—	ON CUSTOMER'S COUNTER
④	IMAGING SYSTEM (UNDER COUNTER)	Ⓢ	110	1,468	17 3/8*	32 1/2*	27*	*INCLUDES 5 1/8" BEHIND, 4" ABOVE AND ALONG CONTAINER.
⑤	POLYDOROS F80 (80 kW) GENERATOR CABINET	Ⓢ	838	2,048**	31 1/2	17 1/8	86 3/4	**DURING OPERATION, 1,195 IN STANDBY MODE
⑥	LUMINOS AGILE (+90/-90) TABLE	Ⓢ	3,793	2,900***	83	88	77	***DURING OPERATION, 512 IN STANDBY MODE
⑦	DCS-1 IC-DVI FLAT DISPLAY & TOUCH DISPLAY ON CEILING SUSPENSION	Ⓢ	419	512	167 1/8	27 7/8	*A	*A - 63" MIN. AND 102" MAX.
⑧	YSIO WALL STAND WITH MOBILE DETECTOR (RIGHT LOADING)	Ⓢ	551	819	30	37*B	83	*B - MAX. IN HORIZONTAL POSITION
⑨	( * ) DOCKING STATION (WALL MOUNTED)	Ⓢ	40	256	20 3/8	7 13/16	16 1/4	WITHIN 11.5 FT. OF WALL STAND
⑩	CEILING RAILS FOR FULLY SYNCHRONIZED TUBE STAND	Ⓢ	59	—	167 3/8	3	3 1/2	SIZE AND WEIGHT PER RAIL
⑪	3M TRANSVERSE BRIDGE AND X-RAY TUBE STAND	Ⓢ	772	853	119 1/4	39	4	
⑫	GRID HOLDER (WALL MOUNTED)	Ⓢ	63	—	25 1/16	6 7/16	21 11/16	SUGGESTED LOCATION
⑬	(MAX) DVI INTERFACE FOR EXAM ROOM DISPLAY (OPTION)	Ⓢ	—	—	—	—	—	CUSTOMER PROVIDES MONITOR
⑭	(MAX) CHARGING STATION FOR MAX DETECTORS (OPTION)	Ⓢ	—	—	—	—	—	ON CUSTOMER'S COUNTER
⑮	(MAX) ACCESS POINT (TOP OF GENERATOR)	Ⓢ	4	—	8	5.5*A	3*B	*A 12 1/4, *B 9" INCLUDING ANTENNAE
Ⓢ	( * ) - THIS ITEM IS ELIMINATED FOR ALL MAX SYSTEMS (MAX) - ITEM ADDED FOR MAX SYSTEMS ONLY		—	—	—	—	—	

# LUMINOS AGILE / AGILE MAX SYSTEM SPECIFICATIONS

## POLYDOROS F80 80kW

### X-RAY GENERATOR POWER REQUIREMENTS

INCOMING POWER:	480 VOLTS, 3 PHASE, 60Hz
CIRCUIT BREAKER:	80 AMPS.
GENERATOR OUTPUT:	80 kW
ALLOWABLE IMPEDANCE:	≤ 0.16 Ω
MAXIMUM MOMENTARY LOAD:	126 kVA
LINE VOLTAGE VARIATION:	± 10% MAX.
PHASE IMBALANCE:	± 2%
FREQUENCY VARIATION:	± 1 Hz

#### NOTE:

ALL INCOMING POWER SUPPLIES, FOR THE SIEMENS EQUIPMENT, ARE TO BE DEDICATED (BACK TO SOURCE) ISOLATED AND INSULATED FROM ANY OTHER EQUIPMENT, SUCH AS, ELEVATORS, GENERATORS, HVAC SYSTEMS, ETC.

A NEUTRAL CONDUCTOR, IF PRESENT, IS NOT USED FOR THE LINE VOLTAGE CONNECTION TO THE SIEMENS EQUIPMENT. IF THE NEUTRAL CONDUCTOR IS PROVIDED, IT SHOULD NOT BE ELECTRICALLY CONNECTED AT ANY POINT IN THE POWER DISTRIBUTION TO THE SIEMENS EQUIPMENT UNLESS SPECIFICALLY REQUIRED. UNINTENTIONAL NEUTRAL TO GROUND BONDS MAY VIOLATE LOCAL AND NATIONAL ELECTRICAL CODES, AS WELL AS CREATE GROUNDING PROBLEMS.

IF AN ON-SITE TRANSFORMER IS REQUIRED TO OBTAIN XP MODALITY OPERATING VOLTAGE, IT MUST BE OF SUFFICIENT CAPACITY AND CHARACTERISTICS TO MAINTAIN SUPPLY VOLTAGE AND IMPEDANCE REQUIREMENTS (TRANSFORMER & CONDUCTORS).

#### ATTENTION:

SIEMENS MEDICAL SYSTEMS, INC. RECOMMENDS THAT THE INCOMING POWER LINES BE ANALYZED WITH RESPECT TO TRANSIENT SURGES AND IMPULSES, SAGS, AND OVERVOLTAGES.

REV 2

## CEILING STAND TECHNICAL DATA

### TRANSPORTING INFORMATION

		SIZE	WEIGHT
TRANSVERSE BRIDGE	3 M	126"L x 32"W x 10"H	419#
	4 M	174"L x 32"W x 10"H	512#
LONGITUDINAL RAILS	4 M	167"L x 3"W x 4"H	59# EACH
	5 M	197"L x 3"W x 4"H	82# EACH
DCS-1/2	4 M	167"L x 3"W x 3"H	441#
X-RAY TUBE STAND (FULLY SYNCHRONIZED)		67"L x 41"W x 53"H	827#
WALL STAND WITH MOBILE DETECTOR (WITH PACKING)		35"L x 93"W x 41"H	898#
MINIMUM DOOR OPENING:		3'-5 3/8"	
MINIMUM CORRIDOR WIDTH:		6'-11"	

### ENVIRONMENTAL CONDITIONS

	IN OPERATION	TRANSPORT
PERMISSIBLE AMBIENT TEMPERATURE (WITH WIRELESS DETECTOR)	59°F TO 82°F	-4°F TO 131°F
PERMISSIBLE RELATIVE HUMIDITY	20% TO 75%	5% TO 95% REV 3

## TABLE TECHNICAL DATA

ENVIRONMENT:	59°F – 95°F OPERATING ROOM TEMPERATURE (NOTE: MAXIMUM 82°F WITH WIRELESS DETECTOR) 20% – 75% PERMISSIBLE RELATIVE AIR HUMIDITY (NON-CONDENSING)
TRANSPORTING/ RIGGING:	<u>UNIT BASE</u> WEIGHT: 1,901 LBS. WITH PACKING 1,610 LBS. WITHOUT PACKING (WITH TRANSPORT CARRIAGE) DIMENSIONS WITH PACKING: 53"L X 36"W X 67"H  <u>UNIT SUPPORT</u> WEIGHT: 2,095 LBS. WITH PACKING 1,517 LBS. WITHOUT PACKING (WITH TRANSPORT CARRIAGE) DIMENSIONS WITH PACKING: 97"L X 52"W X 56"H DIMENSIONS WITHOUT PACKING (WITH TRANSPORT CARRIAGE): 116"L X 32"W X 60"H (WHEELS OUTSIDE) 95"L X 32"W X 60"H (ALL WHEELS INSIDE) 107"L X 32"W X 60"H (TWO WHEELS INSIDE) MINIMUM DOOR OPENING FOR TRANSPORT CARRIAGE: 34"W WITH MINIMUM 8'-3" CORRIDOR WIDTH 48"W WITH MINIMUM 6'-1" CORRIDOR WIDTH

## WIRELESS DETECTOR CONNECTION

OPERATION OF THE WIRELESS DETECTOR CAN BE AFFECTED BY OTHER WLAN DEVICES IN THE VICINITY OF THIS INSTALLATION. TO AVOID ANY CONFLICTS, THE CUSTOMER MUST PROVIDE A LIST OF EXISTING WLAN CHANNELS (FREQUENCIES) OR THE SPECIFIC CHANNEL (FREQUENCY) THEY DESIRE TO BE USED FOR THE WIRELESS DETECTOR.

THE WIRELESS CONNECTION IS ENCRYPTED (WPA2) AND IS BASED ON TWO WLAN STANDARDS, WITHIN WHICH SEVERAL CHANNELS (FREQUENCIES) ARE AVAILABLE:

- 1) 11G STANDARD – OPERATES AT 2.5 GHz
- 2) 11A STANDARD – OPERATES AT 5 AND 6 GHz

THE STANDARD (11G OR 11A) CAN BE SET BY SIEMENS SERVICE VIA THE SERVICE SOFTWARE INSTALLED ON THE IMAGING SYSTEM.

THE WIRELESS CONNECTION IS ONLY USED TO TRANSFER DATA BETWEEN SIEMENS EQUIPMENT AND IS NOT USED TO SEND DATA TO THE CUSTOMER'S NETWORK.

REV 0

FOR REFERENCE ONLY,  
NOT FOR CONSTRUCTION.

# LUMINOS AGILE / AGILE MAX SYSTEM SPECIFICATIONS

## MAXIMUM CABLE DISTANCES BETWEEN COMPONENTS

	CONTROL ROOM MODULE/EQUIP.	IMAGING SYSTEM	GENERATOR	LUMINOS AGILE TABLE	CEILING MONITOR SUSPENSION	CEILING TUBE STAND	DETECTOR WALL STAND	DOCKING STATION
GENERATOR	59'-0"	59'-0"	—	22'-0"	—	32'-0"	36'-0" 52'-6" (NOTE 1)	36'-0"
IMAGING SYSTEM	11'-0"	—	59'-0"	—	59'-0"	—	—	—
LUMINOS AGILE TABLE	—	—	22'-0"	—	—	—	—	—

NOTE 1 — MAX SYSTEM OPTION: 16 METER CABLE EXTENSION MUST BE SELECTED.

THE DISTANCES LISTED ABOVE ARE CALCULATED AS THE MAXIMUM CABLE LENGTH BETWEEN CABLE ENTRY POINTS. DEPENDING ON THE COMPONENT, THE CABLE ENTRY POINT MAY BE IN FLOOR, WALL OR CEILING. VARIOUS ARRANGEMENTS OF COMPONENTS ARE POSSIBLE AS LONG AS THE DISTANCES SHOWN ARE MAINTAINED AND THE SYSTEM FUNCTIONALITY IS NOT ADVERSELY AFFECTED.

## ROOM HEIGHT REQUIREMENTS

### ROOM HEIGHT WITHOUT CEILING STAND

MINIMUM ROOM HEIGHT	RESTRICTION WITH +/-90° ROTATION
8'-2 1/2" TO LESS THAN 10'-4"	<ul style="list-style-type: none"> <li>— MAXIMUM TABLETOP EXTENSION OF 11 1/2" WITH A ROOM HEIGHT OF 8'-2 1/2"</li> <li>— DIGITAL IMAGING TOWER CANNOT BE IN THE EXTREME HEAD-END POSITION WITH A ROOM HEIGHT OF 8'-2 1/2"</li> <li>— RESTRICTIONS IN THE MOVEMENT RANGE WITH AUTOMATIC CORRECTION OF MOVEMENT.</li> </ul>
10'-4" OR GREATER	NO RESTRICTIONS

### ROOM HEIGHT WITH CEILING STAND

	USABLE TABLE HEIGHT AT 45° SID (3)	
MINIMUM ROOM HEIGHT	8'-9 1/8"	2'-6"
MINIMUM ROOM HEIGHT FOR 60° SID TO TABLE (3)	9'-7" (1)	3'-4"
MAXIMUM ROOM HEIGHT WITHOUT THE TUBE STAND TELESCOPE EXTENSION	9'-5 5/8" (1)	3'-2"
	9'-10 7/8" (2)	3'-7 3/8"
MAXIMUM ROOM HEIGHT WITH THE TUBE STAND TELESCOPE EXTENSION	10'-1 1/2" (1)	3'-2"
	10'-7" (2)	3'-7 3/8"

1) UPRIGHT (0°) EXPOSURES ARE POSSIBLE AT LOWEST POSITION OF WALL STAND.

2) UPRIGHT (0°) EXPOSURES ARE NOT POSSIBLE AT LOWEST POSITION OF WALL STAND.

3) WITH SECOND PLANE OF CEILING STAND.

## REMOTE SYSTEM DIAGNOSTICS

SIEMENS REMOTE SERVICES (SRS) REQUIRES A CONNECTION BETWEEN THE SRS REMOTE SERVER AND SIEMENS SYSTEMS VIA REMOTE LOCAL AREA NETWORK ACCESS, TO ENSURE THE UPTIME OF YOUR SYSTEM. A CUSTOMER VPN CAPABLE FIREWALL OR OTHER VPN APPLIANCE IS PREFERRED.

## FOR MORE INFORMATION

FOR MORE DETAILED PLANNING REQUIREMENTS FOR THIS SYSTEM, SEE THE TYPICAL FINAL DRAWING SET NUMBER: 11028